#whdentalwerk (f) (☉) (m) (►) video.wh.com Now at your dealer or **wh.com**






Incredible inside & outside

New Lara sterilizers Ready for the future



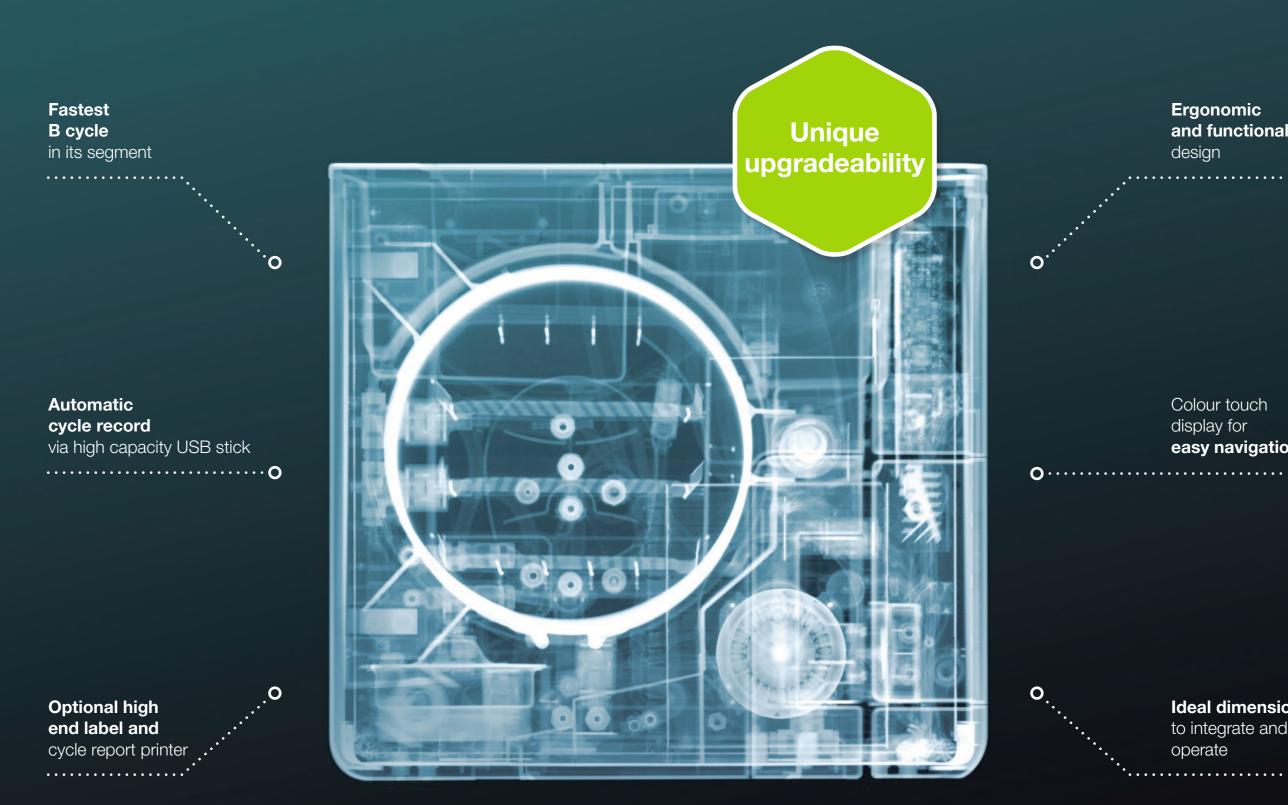

X-ray image by Nick Veasey for W&H. Find out more about the incredible X-ray-work by artist Nick Veasey at nickveasey.com

Ergonomic and functional design

Colour touch display for easy navigation

Ideal dimensions

to integrate and operate •

Discover simple operation & reliable traceability


The Lara colour touch screen has been selected and designed to make your daily work easier, faster and more efficient. Lara also offers a comprehensive traceability system customized to your needs.

Easy to use

The intelligent menu structure of the new Lara sterilizers ensures extraordinarily high operating comfort: you can browse to the preset programmes quickly and intuitively using the large colour touchscreen.

High capacity data logger

A high capacity USB drive automatically records the cycle reports throughout Lara's entire service life. Optionally available: the Lara label and cycle report printer with which you can document the traceability without an additional computer or separate software.

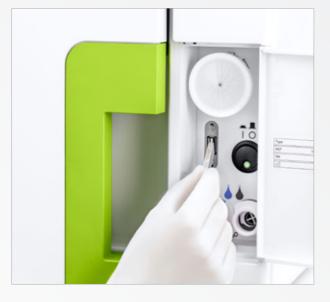





Easy filling A high capacity integrated funnel prevents water splashes during the filling process.



Effortless tank access The tank cover can be easily removed without tools to clean the tank.

lara

(WSH)



User identification

features such as user identification with Activation Codes: the user who has loaded/unloaded the sterilizer can be digitally recorded.





You can also acquire additional Users are identified by a 6-digit PIN.

Load release confirmation

Also available as an upgrade: the approval confirmation. This option digitally records the name of the end user and their confirmation of the successful sterilization and/or test cycle.

Automatic water filling An automatic water filling

valve allows you to connect a demineralization system. Manual filling and draining are no longer required.

Discover unique upgradeability

Whatever occur in the future, with the new Lara sterilizers you're well prepared. Effectively, their innovative upgradability system offers you a cost-efficient, rapid and customized opportunity to activate additional features to reflect your practice needs or comply with future requirements.

The four additional features at a glance:

Even faster cycle times by

saves time and energy.





Discover outstanding performance

The new Lara sterilizers offer you one of the fastest type B in its performance class right from the standard version. And you can improve them even by activating the additional features Eco Dry + and Fast Cycle.

Eco Dry + automatically adapts the drying time to the load and, in doing so, saves time and energy. Additionally available via an Activation Code: a type S fast cycle for unwrapped, semi-critical products.

eco dry+

Time saving

Type B cycles process and perfectly dry 2 KG OF LOAD IN LESS THAN 40 MINUTES.

Extends service life

Automatic adjustment of the drying time depending on the load reduces the heat exposure time. This extends the service life of your instruments.

Energy saving

The optimization of the drying time means reduced energy consumption. Making Lara the »green solution«.



Reprocessing in less than 23 minutes!



Disinfect your instruments before cleaning and servicing by wiping them.

10"

Wipe

disinfection

in 1 minute

Sterilization in 20 minutes With the Lara fast cycle, the instruments can be removed from the sterilizer ready for use after only 20 minutes.



Technical

data

lara



Sterilization cycle		B-Universal 134		B-Prion 134		B-Universal 121	Fast cycle Optional
Temperature:		134°C		134°C		121°C	134°C
Holding time (minutes):		4'		18.5'		20.5'	3.5'
Total cycle duration (minutes) including complete drying:	Lara 17	Automatic from 25' to 51'	37'	Automatic from 40' to 65'	52'	41' to 73'	20° to 22°
	Lara 22	Automatic from 28' to 55'	39'	Automatic from 40' to 68'	54'	42' to 76'	20' to 22'
Lara 17 loads: Lara 22 loads:		maximum 4.5 kg/typical 2 kg/maximum porous 1.5 kg maximum 5.5 kg/typical 2 kg/maximum porous 2 kg			max. unwrapped 2 kg max. unwrapped 2 kg		
Lara 17/22 containers load:		maximum 9 kg					
Test cycles:		Helix/Bowie and Dick/Vacuum					

Items supplied with reversible rack, 3 aluminium trays, drain tube, tray holder, 8 GB USB drive, door opening tool, factory test

CE.... CE....

Lara 17	Lara 22						
Lara RIL-117	Lara RIL-122						
171	22						
200–240 V AC; 50/60 Hz; 10 A	200–240 V AC; 50/60 Hz; 10 A						
2.0–2.4 kW							
465 x 452 x 646 mm							
42.5 kg	44 kg						
4.8/4.8							
7 to 15 cycles							
195 x 195 x 312 mm	195 x 195 x 400 mm						
2 USB ports (5 as option), integrated automatic water filling valve							
Lara sterilizers were designed, certified and validated with the most stringent directives and standards:							
EN 13060 Small steam sterilizer	IEC 61010-1 Safety requirements						
IEC 61326-1 Electromagnetic compatibility	IEC 61010-2-040 Specific requirements for steam sterilizer						
	Lara RIL-117 17 I 200-240 V AC; 50/60 Hz; 10 A 2.0-2 465 x 452 42.5 kg 42.5 kg 42.5 kg 7 to 15 195 x 195 x 312 mm 2 USB ports (5 as option), integr						

IEC 61770 Electric appliances connected to the water mains

> The sterilizer can be validated according to EN ISO 17665-1

Safe, complete reprocessing thanks to W&H

Sterilization, hygiene and maintenance are increasingly important in the dental practice. W&H offers best practice hygiene protocols as a complete solution for processing your straight and contra-angle handpieces and turbines.



W&H offers a wide range of accessories for optimization of the sterilization process. Minor extras lead to major results. Work in the practice is easier and your team saves time.

disinfectors

Maintenance devices

Sealing devices

It is all about **your demands**

For years W&H has been known as pioneer in terms of sterilization. Our biggest aim is to develop innovative, reliable and fast sterilizers, that fit perfectly into the individual practice environment and workflow. With our range of sterilizers we are able to offer the perfect sterilizer for each demand.

Find out more about our sterilizer family: now at your local dental dealer or at **wh.com**



W&H Sterilization Srl. Via Bolgara, 2 24060 Brusaporto (BG), Italy t +39 035-66 63 000 f +39 035-50 96 988 office.sterilization@wh.com wh.com Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.